



The 10 LPM "Hi-Flow Pro."
The low-cost modality for high-flow patients.

THE FIRST AND BEST 10 LPM CONCENTRATOR JUST GOT BETTER

E-Z TO READ FLOW INDICATOR

E-Z TO CLEAN EXTERIOR

E-Z TO SERVICE




A Reliable, Cost-effective Way To Serve The High-flow Patient

Patented ATF® module (Advanced Technology Fractionator) eliminates more than 200 troublesome parts found in conventional oxygen concentrators!


- No failure-prone solenoid valves
- 40 fewer pneumatic connections
- 35 fewer electrical connections
- No electronics for valve sequencing
- No pressure-reducing regulator

One mechanically integrated module, with a single moving part made of space-age composites

- Maintenance-free
- Self-cleaning
- Invulnerable to wear

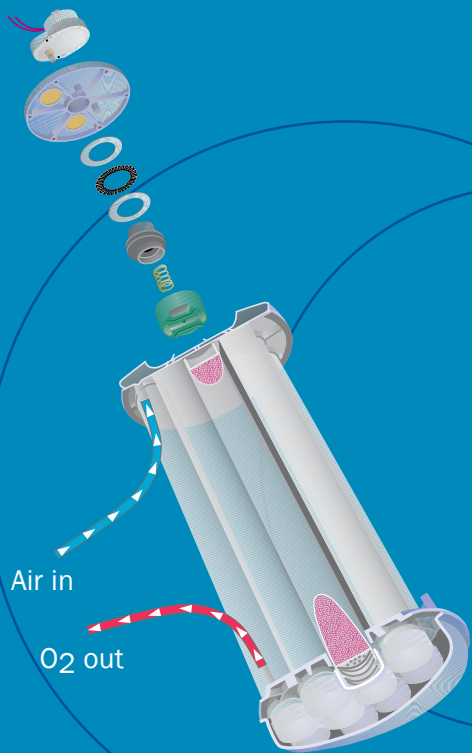

A compressor built for steady-state operation is matched to the ATF's constant low-pressure operation. The result? Elimination of the pressure swings that can lead to costly compressor rebuilds.

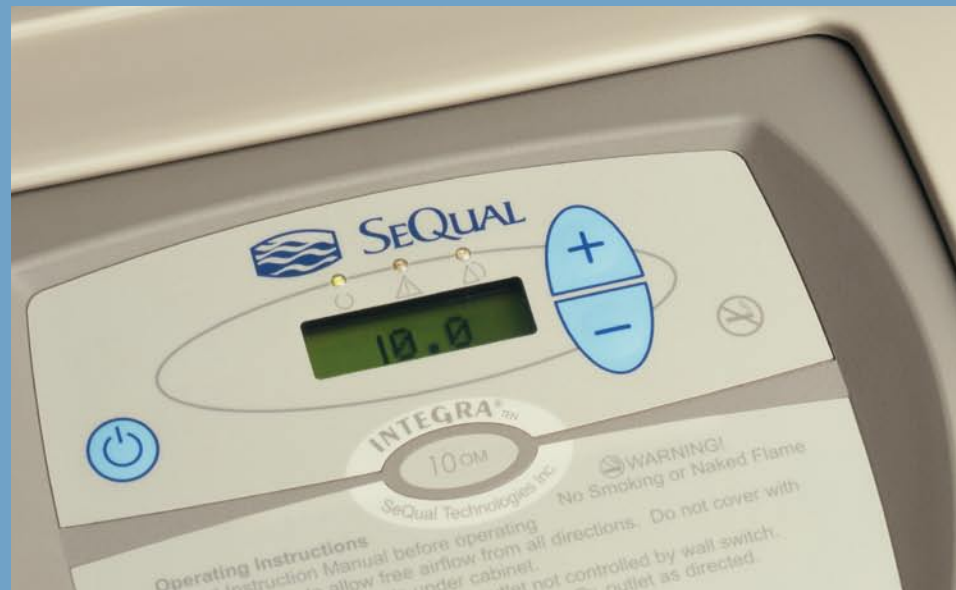

WHY SETTLE FOR ANYTHING LESS THAN A PERFECT 10?

Integra Ten is the "High-Flow-Pro" with over 5 years of proven reliability. No conventional concentrator can match SeQual's patented ATF® performance. The cool running Integra is the perfect choice for providers and users, and now easier than ever to operate, clean, and service.

For the Dealer/Provider . . .

- The most reliable, cost-effective system for serving high-flow patients.
- Eliminates the costly oxygen deliveries needed with alternate modalities.
- Easy-to-read flow indicator eliminates patient concerns and phone calls to confirm flow rate.
- Patient cannot damage unit by setting flow beyond maximum capacity.
- Internal 9 volt battery, with provider accessible reserve gauge, powers alarm if electric service is interrupted.
- An hour meter can be easily accessed and viewed by service technician.
- Essentially the same size, weight, and power consumption as a 5 LPM machine.

For the User . . .

- Simplist keypad controls: On/Off and O₂ increase/decrease arrows.
- Accurate, easy-read LCD flow indicator in ½ L increments.
- New flat panel cleans with a quick wipe.
- Alerts user if O₂ hose becomes restricted.
- Delivers up to 10 LPM of continuous O₂ with the security and convenience only a concentrator can provide.
- Single, easy-to-clean cabinet intake filter is highly visible and easy to reach and remove for periodic cleaning.
- Compact cabinet with integrated handle and four swivel casters makes movement easy.

Adjustments have never been simpler than with the large, easy-to-read LCD flow indicator. Flat and easy to wipe clean, it has no buttons or knobs that can get knocked off or damaged. The result? Nuisance and unscheduled service calls have never been fewer.



A "zero-clearance" air intake located atop the Integra™ prevents overheating if the unit is pushed against walls or draperies. An ergonomically designed handle has a sure grip with hand-fitting comfort.



A Velcro power cord holder simplifies storage and transport.



A DISS-fitting cabinet design accommodates extra large humidifiers.



Large-diameter steel-forked casters hold up under rough handling and make clearance, mobility easy. Rear wheel brakes ensure stability.



Integra™ Specifications

Performance

Flow rate: Variable 1 to 10 LPM (1/2 liter increments)

O₂ concentration: 91%±3% @ sea level (1-10 LPM flow)

O₂ pressure: 7 PSIG (48 Kpa) nominal

Patient Controls

Control panel-mounted On/Off switch

Flow adjustment buttons

Alarms

Power failure

O₂ concentration (optional)

Irregular pressure

General

Size: 26.0 H x 14.7 W x 19.1 D (inches)

58.4 H x 39.4 W x 47.0 D (cm)

Weight: 57 lbs. (26 Kg)

Approvals

Class II Equipment

Double Insulated ■

Type BF Applied Part □

CE Mark (Model 6323A-OM-10 only)



200972



Customer Service Department

11436 Sorrento Valley Road

San Diego, CA 92121-1306

Telephone: 800.826.4610

858.202.3100

Toll-Free Order Fax: 877.786.9943

E-mail: medicalsales@sequal.com

Web Site: www.sequal.com

Model Selection Guide

Part No. 9622 Integra E-Z Model 6323A-10

Without O₂ monitor

110 volt; 60 Hz

English

Part No. 9620 Integra E-Z Model 6323A-OM-10

With O₂ monitor

110 volt; 60 Hz

English

Part No. 9620-1 Integra E-Z Model 6323A-OM-10

With O₂ monitor

110 volt; 60 Hz

English/French

Part No. 9620-3 Integra E-Z Model 6323A-OM-10

With O₂ monitor

110 volt; 60 Hz

Spanish/Greek

Part No. 9623 Integra E-Z Model 6323A-10

Without O₂ monitor

220 volt; 50 Hz

English

Part No. 9621 Integra E-Z Model 6323A-OM-10

With O₂ monitor

220 volt; 50 Hz

English

Part No. 9621-1 Integra E-Z Model 6323A-OM-10

With O₂ monitor

220 volt; 50 Hz

English/French

Part No. 9621-2 Integra E-Z Model 6323A-OM-10

With O₂ monitor

220 volt; 50 Hz

German/Dutch

Part No. 9621-3 Integra E-Z Model 6323A-OM-10

With O₂ monitor

220 volt; 50 Hz

Spanish/Greek

Specifications subject to change without notice.